EPilaser

MODE D'EMPLOI | MANUAL | MANUAL
MANUALE | HANDBUCH | HANDLEIDING

EPILADY®

Français	1
English	24
Español	47
Italiano	69
Deutsch	92
Nederlands	115

TABLE DES MATIÈRES

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TABLE OF CONTENTS

English

PREAMBLE	25
INTENDED USE OF THE EPILASER (INDICATIONS)	26
WHO SHOULD NOT USE THE EPILASER? (CONTRAINDICATIONS)	26
Warnings	27
Important Safety measures	29
1. Welcome	30
1-1 Introduction	30
1-2 Features	31
1-3 Box Contents	31
2. Device Description	32
2-1 Overview	32
2-2 Charging the device	32
2-3 Switching the device on/off	34
3. Getting Started	34
3-1 Preparation	34
3-1-1 Check your skin tone and natural hair color	35
3-1-2 Prepare your skin	36
3-1-3 Skin test on small area	37
3-2 Treating your hairs	37
4. After the treatment	41
4-1 Cleaning the device	41
4-2 User Tips	41
Troubleshooting	44
Specifications	45
Notices	46

PREAMBLE

Important: Please read this manual before using this device.

Users are required to familiarize themselves with all the safety information and instructions before using this device. This manual should be saved for future reference and use.

Description of device markings.

5V  1000mA : Power

Laser class 1C: energy source via 3 laser diodes

 **Warning :**

The laser output can be dangerous if the user does not use the device in accordance with the instructions.

NOTE: The device has been tested according to a protocol defined on page 35 of this manual. More frequent use is not recommended and may be a risk.

 The device is equipped with a battery that is not user-replaceable and is intended to be recharged by a removable power supply unit.

 **Warning :** This device contains a battery that can only be replaced by qualified personnel.

WARNING : This device contains a battery that can only be replaced by qualified personnel.

 : Product complies with European requirements in accordance with the legislation in force to date.

 : This device cannot be disposed of automatically, please contact the sorting center closest to your home.

INTENDED USE OF THE EPILASER (INDICATIONS)

Epilaser is an over-the-counter device designed to be used in conjunction with shaving for long-lasting hair removal through periodic treatments. The Epilaser is also designed for the permanent reduction of hair regrowth, defined as a stable, long-term reduction in the number of hairs after a treatment period. Permanent hair reduction is defined as a stable, long-term reduction in the number of hairs growing back, counting 6, 9, and 12 months after the end of the treatment period.

It is designed for use only on light to dark skin (see “Check your skin tone and hair color” on page 35) and on naturally dark hair. The Epilaser will work, but will not be effective on white, grey, blond, or red hairs.

WHO SHOULD NOT USE THE EPILASER? (CONTRAINDICATIONS)

Skin colors not compatible: Type V and Type VI

DO NOT USE the Epilaser if your skin tone corresponds to one of the colors marked “not compatible” on the skin tone safety chart (see “Check your skin tone and hair color” on page 35). If you use the Epilaser on dark, very dark, or black skin, indicated as “not compatible” for treatment, the skin may absorb too much laser light causing burns, discoloration, or scars.

NOTE: The Epilaser is equipped with a skin sensor. If the sensor determines that its use is “not compatible” with your skin tone, the treatment will not take place.

DO NOT USE the Epilaser if you are under the age of 18. Do not use the Epilaser on yourself or on anyone under the age of 18. Using the Epilaser on yourself or by children under the age of 18 has not been studied and could result in serious injury.

DO NOT USE the Epilaser if you have a light hair color: white, red, or similar to blond, white, red, or gray. For the laser to be effective, you must have naturally dark hair in the area to be treated. Light-colored hair pigments won't absorb enough laser energy to destroy the hair follicle.

Warnings

1. DO NOT USE the Epilaser on the following areas of your body:

- On or around the eyes, eyebrows, or eyelashes. This may cause permanent eye damage. Anywhere on or around the forehead, temples, nostrils, ears, nipples, genitals, or around the anus. These areas may have more sensitive skin, a darker complexion and/or greater hair density. Use of the Epilaser in these areas has not been tested, and could cause serious skin damage.
- On parts of your skin that are darker than your usual skin tone, such as tattoos, brown or black spots (e.g., large freckles, birthmarks, or moles), tanned skin, or red or swollen skin. Dark skin can absorb too much laser light, which can burn the skin.
- On damaged skin, including open wounds, cold sores, sunburn, or areas where you've had a peeling or other skin resurfacing procedure within the last 6-8 weeks (check with your doctor). These conditions can make the skin more sensitive to the laser, resulting in burns.
- On any area where you'd like your hairs one day to grow back.

DO NOT USE the Epilaser if your skin tone is identified as "not compatible" (See "Check your skin tone and hair color" on page 35).

If in doubt, please see "Check your skin tone" on page 35. Darker skin tones can absorb too much light, which can damage the skin.

2. DO NOT USE the Epilaser if you are pregnant or breastfeeding.

The Epilaser has not been tested on pregnant women. Hormonal changes may increase skin sensitivity and the risk of skin lesions.

3. DO NOT point the laser at your eyes.

The laser could cause permanent and irreversible damage

4. DO NOT USE the Epilaser for any purpose other than hair removal.

5. DO NOT USE the Epilaser if it seems to be damaged. DO NOT attempt to open or repair the Epilaser.

These actions increase the risk of skin or eye injury. Please contact your dealer's customer service department for assistance if your Epilaser seems to be damaged.

6. Always follow the basic safety precautions for electrical appliances:

- DO NOT get the Epilaser, power cord, or adapter wet.
- DO NOT use the Epilaser while bathing, and do not place it near a bathtub or sink.
- DO NOT USE a power adapter with a damaged plug or cord.
- ALWAYS keep the cord away from heated surfaces. Failure to observe these precautions could result in electrocution or electric shock.

7. DO NOT drop the Epilaser or any of its components, as this may damage the unit.

Important Safety Measures

- 1.** Humidity: 25%–80%, Storage temperature: –20°C to 60°C(–4F to 140F), Operating temperature: 0–50°C(32F to 122F).
- 2.** DO NOT insert anything into the treatment window opening or any other part of the laser.
- 3.** DO NOT get the Epilaser and/or power cord wet. If they get wet, do not use the Epilaser or the power cord. Performing a treatment with a wet Epilaser (or component) can cause electrocution or electric shock. If any part accidentally gets wet, dry it thoroughly before using it again.
- 4.** Use only original or certified cords. DO NOT use the power cord with a damaged plug or cable.
- 5.** Dropping the Epilaser or any of its components may damage the device.
- 6.** Keep the storage temperature between 0° C to 35° C.
- 7.** Keep the device away from excessive heat, such as sunshine, fire, etc. Likewise, keep it away from excessively cold temperatures.
- 8.** Do not use the device outdoors during thunderstorms.
- 9.** Keep out of reach of children. It contains small parts that may pose a choking hazard.
- 10.** To avoid damage to parts or internal circuitry, do not use your device in dusty, smoky, damp, or dirty environments.
- 11.** Do not use, store, or move the device near explosives or flammable objects. Using your device in these environments increases the risk of explosion or fire.
- 12.** Clean the unit with an alcohol wipe. Do not use water or corrosive products/oil to clean it.

CHAPTER 1 WELCOME


1-1 Introduction

The Epilaser is a permanent hair reduction laser that gives you long-lasting results in the comfort of your own home, unlike temporary solutions like waxing. The Epilaser works in the same way as professional permanent laser hair reduction. The laser beam heats the dark pigment inside the hair and deactivates the hair follicle, giving you permanent hair reduction after several sessions. It is only effective on natural dark hairs (if in doubt, refer to page 35).

The built-in camera and sensors enable the Epilaser to automatically find and treat all hairs in the area to be treated in a matter of seconds. Simply hold the Epilaser over the area to be treated and move on to the next area once the treatment is completed.

Where to use the laser?

The Epilaser can be used on the following areas of the body: Women's facial hairs (cheeks and below, such as chin and jawline), arms, armpits, hands, legs, feet, bikini line, and back.

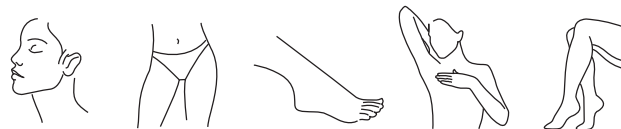Warning

It is dangerous to use the Epilaser on the following areas of the body: Eye area (eyes, eyebrows, and eyelashes), forehead, ears, nipples, genitals, and anus.

1-2 Features

1. Permanent results: Permanent hair reduction and smooth skin.


2. Easy to use: To use the Epilaser, simply place it on the skin and press a button.

3. Safe and efficient: Uses professional laser technology, delivering the ideal amount of light with each pulse for safe and efficient use.

4. Practically pain-free: the Epilaser provides a gentle, practically pain-free treatment.

5. Convenient: Avoid appointments and treat your skin in the comfort and privacy of your own home.

1-3 Box Contents

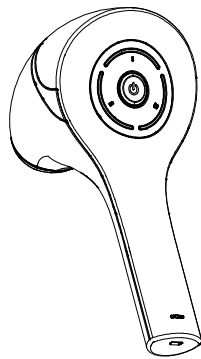
EPILASER

USB cable


AC power adapter

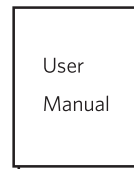User
Manual

CHAPTER 2 DESCRIPTION OF THE DEVICE


2-1 Overview

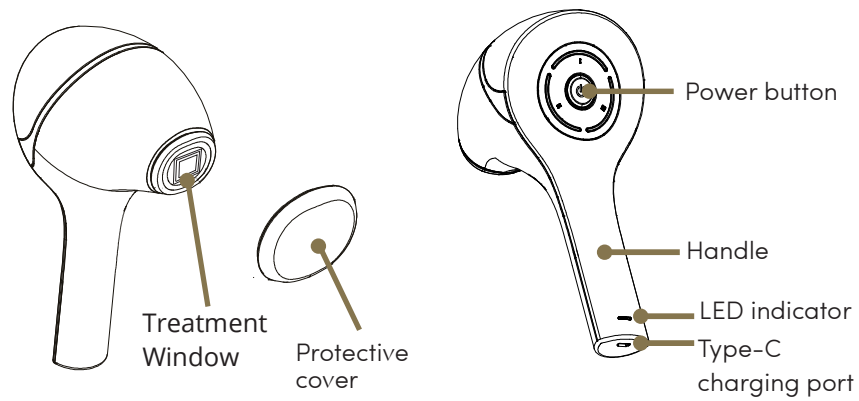2-2 Charging the Device

Notes


- Before charging your device, make sure the charging port is clean and dry.
- Using other brands of chargers may result in slow charging and/or overheating. Use the adapter and USB charging cable supplied.
- To charge your Epilaser, connect one end of the supplied power cord to the Epilaser's charging port and the other end of the cable to a USB power source.

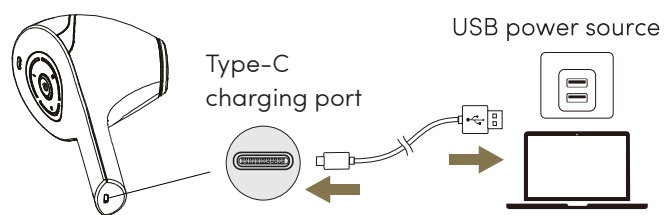CHARGING IN PROGRESS
 ● Steady red

FULLY CHARGED
 ● Steady green

Note: Use only the adapter and USB cable supplied with your Epilaser.

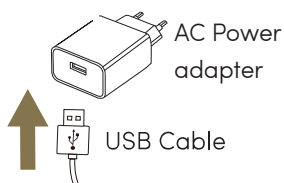
For the Epilaser sold in the European region, first connect the USB end of the included USB charging cable to the power adapter to create a charger. Then connect the charger to the charging port of the Epilaser and to a wall outlet.



- It takes about 3 hours to fully charge the Epilaser.
- When the Epilaser has finished charging, disconnect the power cord and USB power source.
- If you don't want to use your Epilaser for a long time, we recommend re-charging every 3 months to preserve battery life.

2-3 Switching on/off

1 Remove the protective cap from the Epilaser.

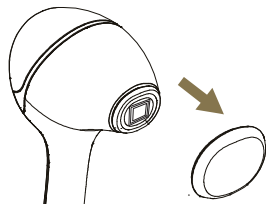2 Press and hold the activation button for 2 seconds to turn the Epilaser on or off.

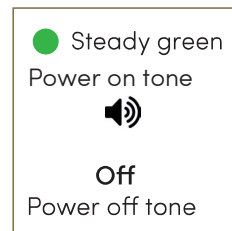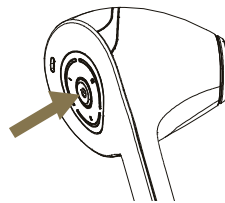CHAPTER 3 GETTING STARTED

3-1 Preparation


Follow these instructions to prepare your treatment.

1. Do not perform a peeling or any other skin resurfacing procedure six weeks or less before your first treatment.
2. Stay out of the sun to avoid sun exposure before and after treatment. When you go out, apply sunscreen with a high protection factor (SPF30).
3. Avoid self-tanning creams that darken your skin.
4. Avoid other hair removal methods. Hair removal with tweezers, waxing, and electrolysis can disrupt the hair follicle and should be avoided for at least four weeks before the treatment.

3-1-1. Check your skin tone and natural hair color

- Check your skin tone

To determine whether this laser is safe to use, compare the skin color of the area(s) you wish to treat with the skin colors shown below. The Epilaser is intended for use only on skin tones from Type I (white) to IV (light brown).

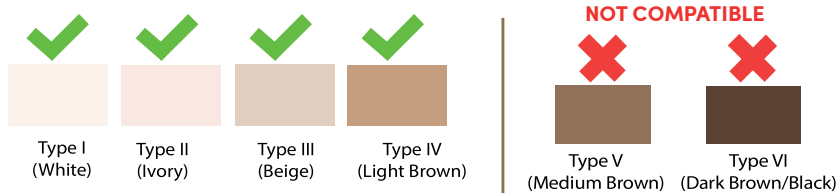
Warning

DO NOT USE the device if the skin tone of the area to be treated corresponds to Type V (medium brown) or Type VI (dark brown/black). Dark skin can absorb too much laser light, which can damage the skin.

- Check your natural hair color

Epilaser works best on naturally dark hairs (black, dark brown, medium brown, or light brown). The Epilaser will work, but will not be effective on white, grey, blond, or red hairs.

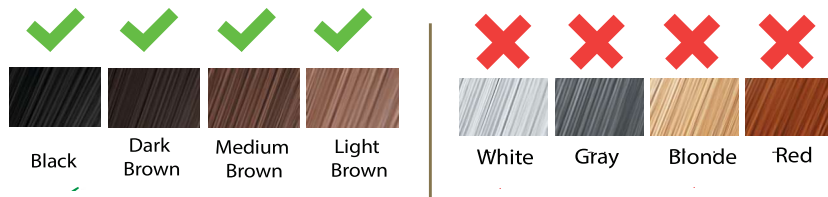
3-1-2 Prepare your skin

STEP 1. Trimming and shaving

We recommend trimming and shaving the hairs the day before the laser treatment. Trimming or shaving hairs above the skin to avoid superficial skin damage from burnt hairs, but leave the hair shafts intact below the surface.

The hairs must be visible at the time of treatment (from the shortest to the longest 2 mm hairs). Do not treat hairs longer than 2 mm. Any hair that is too long or too thick will burn.

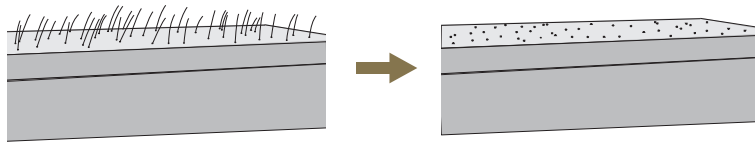STEP 2. Clean the treatment area

Remove all cosmetics, lotions, and creams that could block the laser light. It may be preferable to use the laser after a shower or bath to ensure that the skin is clean.

STEP 3. Make sure that the treatment area is completely dry before starting.

3-1-3 Skin test on a small area

1. Prepare your skin.
2. Turn on the Epilaser
3. Follow the instructions in the “Treating your hairs” section on pages 37 to 40 to treat a single area.
4. Move the laser to the next treatment area and test it with a second treatment area.
5. Please test up to 3 zones.
6. Wait 24 hours and check the tested areas. If there’s no serious discomfort or side effects regarding your skin, the test is successful and you can start treating your hairs.

If you are unsure of the test result, do not use the Epilaser and contact the vendor’s customer support service before continuing to use it.

3-2 Treating your hairs

Note: The laser beam impacts the follicle (“hair root”) beneath the skin’s surface. On impact, the beam heats the follicle, weakening it and preventing future growth. To be effective, each hair follicle must be treated during its active growth period. As hair grows in cycles, several treatments are required. The hair will fall out on its own after a while, or you can remove it once the treatments have been completed.

Tip: For optimum results, we recommend using your Epilaser at maximum level (level 3), with 2 flashes per treatment and 2 sessions per week for 12 weeks. If the maximum level proves too sensitive for the area treated, reduce to level 2 or even 1, while maintaining 2 flashes/treatment and 2 treatments per week.

PLEASE NOTE THAT THE REDUCTION IN POWER MAY EXTEND THE DURATION OF THE TREATMENT.

 **Warning**

Caution - Use of controls or adjustments or performance of procedures other than those specified herein may result in hazardous radiation exposure.

STEP 1. Turn on your Epilaser

Remove the protective cover from the Epilaser and turn it on. Repeatedly press the Power button to select your desired treating level. For optimum results, we recommend using your Epilaser at High level (level 3) at first. Once the laser is turned on, your device is ready for use.

NOTE: If you do not use the device for 2 minutes, the Epilaser will turn off automatically.

STEP 2. Target and treat hairs


Keep the flush treatment window on the area of your skin that you wish to treat. If the skin tone sensor detects that the skin type is not suitable for treatment, the laser will not perform the treatment. Simply hold the device in place on the skin's surface, and it automatically applies the laser treatment. The integrated camera and sensors find and target hair follicles. The laser treats hairs found in the target area.

When the laser is processing, the LED indicator blinks green.

When all the hairs in the target area have been treated, a sound signal indicates that you can remove the device from your skin and move on to the next treatment area.

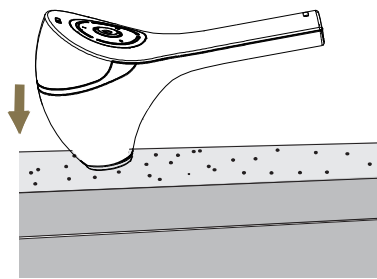If the Epilaser is moved too soon, the indicator light blinks red and you will hear a sound signal, indicating that the treatment has not been completed.

If this happens, simply reposition the device, hold it against the skin and don't move it until it notifies you that it has completed the treatment.

The hairs are treated		Blinking green
The treatment is completed		A sound signal to indicate the end of the treatment
The treatment is not completed		Blinking red
		Failure sound signal

Please note the following points when treating unwanted hairs:

The Epilaser emits several pulses from the treatment window when the area is treated. Once the treatment is complete, the Epilaser must be removed from the skin and repositioned. The Epilaser will only apply one series of pulses if it is not removed from the skin.

Treatment tips

Proper use of the Epilaser requires some practice. At first, it may take some time to hold the device correctly so that it doesn't lose contact with the skin. It may be helpful to use both hands, or to rest your arm on a flat surface for greater stability (It is not necessary to press hard against the skin). Once you have more experience with the Epilaser, treatments will be easier - and you can even have fun!

To help you position the Epilaser correctly on your skin, you can use a mirror. Sitting in front of the mirror will enable you to target the desired treatment area before applying the laser to your skin.

STEP 3. Reposition the device and repeat on the next area

After processing an area, move and reposition the processing window on another area. Follow step 2 "Targeting and treating hairs" (see page 38) for each area you wish to treat.

STEP 4. Turning the device off

When the treatment is completed, turn the device off and clean the Epilaser thoroughly after each use (see "Cleaning the Epilaser" on page 41). Replace the protective cover over the treatment window. Clean the treatment area and, to reduce any discomfort, you can put ice to it.

When used following the instructions, the Epilaser is an over-the-counter device for personal use that provides permanent hair reduction. Results may vary from person to person, depending on factors such as the area of the body treated and compliance with the treatment instructions.

Important:

For best results, treat continuously for 12 weeks, 2 treatments a week. We recommend a double flash/treatment, i.e. 4 flashes per week.

Depending on your body hair and the laser intensity used, the treatment may be extended beyond 12 weeks.

CHAPTER 4 AFTER THE TREATMENT

4-1 Cleaning the Epilaser

You must keep your Epilaser clean and dry, especially the treatment window area. If cleaning is necessary, turn the laser off. Clean the outside of the device with alcohol wipes. It is also advisable to regularly remove dust from its gaps.

- In particular, clean the areas that have touched your skin, especially around the treatment window.
- Take care not to damage the inside or around the treatment window.

Note: DO NOT moisture the unit when cleaning, as this may cause electric shock or electric shock hazard.

4-2 User tips

After treatment, you may have a feeling that varies depending on the person and the treatment area. You may feel nothing, or have a warm, tingling feeling. Some users may see a slight redness, which generally disappears within 24 hours.

We recommend avoiding exposure to artificial light sources, such as tanning booths, sunbeds, sunlamps, or tanning beds.

For 24 hours after treatment, avoid any activity or product that could irritate your skin. This includes using jacuzzis or saunas, wearing tight-fitting clothes, using bleaching creams, perfumed products, or exfoliating products. The above-mentioned products may cause irritation, such as itching and redness, in the treated areas.

Stay out of the sun to avoid sun exposure before and after treatment. When you go outside, apply sunscreen with a high protection factor (SPF30).

Treatment period by sessions

After your first treatment session, there will be no visible change regarding the hairs. Laser hair removal does not cause immediate hair loss. Treated hairs will eventually fall out on their own. Continue with your normal depilatory routine (waxing, shaving, etc.), if you want to.

With further treatments, a gradual reduction in the number of hairs can be seen. Regrown hairs are also lighter and finer.

After the session:

Hair follicles that have been deactivated by the Epilaser will not grow back. Perform occasional touch-ups or additional treatments as needed.

NOTE: Laser hair treatment does not immediately remove treated hairs. As the permanent result, hair follicles will be deactivated on impact and will not grow new hairs.

TROUBLESHOOTING

Problem	Solution
My Epilaser won't turn on.	Check whether the Epilaser battery is low. If so, recharge it.
While charging my Epilaser, the LED indicator doesn't turn on.	Check that your Epilaser is correctly connected. Confirm that you are using the type C power cord supplied with the Epilaser kit, and that the power cord is connected to a working USB source.
My Epilaser was working fine, then it stopped working.	Check whether the Epilaser battery is low. If so, recharge it.
Everything is done correctly, but the Epilaser doesn't work.	Check the following: 1. Check whether the laser has been turned on. If not, press and hold the activation button until the green LED turn on. 2. Check whether the integrated skin tone sensor has detected a skin tone identified as "not compatible". If this is the case, the skin tone in this area is too dark for laser treatment. 3. Check that the treatment window is in full contact with your skin.

Specifications

Specifications are subject to change without notice.

Model Name	Epilaser (permanent hair reduction laser)
Model Number	EP-720-10
Laser source	Diode laser
Laser Wavelength	808 +/- 20% nm
Spot Size	1 mm ²
Energy applied to the skin	24 J/cm ²
Pulse duration	0.1 - 0.2 sec
Size of treated area	18 mm x 14 mm
Sensor	4 X Skin sensor
Safety system	4 Skin sensors+Hardware limitations+Software procesing unit (ARM cortex M7 arc)
Charging port	Type-C port
Power Input (Charging)	5 V/2A
Battery type	3.7V 2600mA rechargeable Lithium-ion Battery
Charging time	Up to 3 hours
Weight	390 g/ 13.75 oz

Notices

This symbol means that the product should not be disposed of with household waste, but must be handed over to an appropriate collection center for recycling.

Comply with local regulations and never dispose of the product and rechargeable batteries with normal household waste. Proper treatment of old products and rechargeable batteries helps avoid negative consequences for the environment and human health.

Disclaimer

All information, design, and specifications contained in this manual were correct at the time of publication. However, as the product is subject to continuous updates and modifications, the final device may differ slightly. The contents of this manual is subject to change without notice, and Epilady will not be liable for errors contained herein or for direct or indirect damages arising from product performance or use of this manual.

Epilady® IS A REGISTERED TRADEMARK OF EpiladyUSA, Inc. and USED UNDER LICENSE BY Epilady 2000 LLC.

©2023 EpiladyUSA, Inc. ALL RIGHTS RESERVED AND ANY USE OF THESE MARKS BY Epilady 2000 LLC. IS DONE UNDER LICENSE.

All rights reserved. Epilaser, Epilady and the Epilady logo are registered trademarks of Epilady 2000 LLC.

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Epilady® IS EEN GEDEPONEERD HANDELSMERK VAN EpiladyUSA, Inc. EN GEBRUIKT ONDER LICENTIE DOOR Epilady 2000 LLC.
TEpilaser, Epilady en het Epilady-logo zijn geregistreerde handelsmerken van Epilady 2000 LLC.